

Bacterial WGS training : Exercise 3


Title	SNP-based bacterial outbreak investigation.
Training dataset:	
Questions:	<ul style="list-style-type: none">• Do I have the needed depth of coverage?• Have I chosen the correct reference?• How do I create a SNP matrix? How many SNPs do I have?• How can I visualize my phylogenetic tree? Which problems can I encounter?• Which strains belong to the outbreak?
Objectives:	<ul style="list-style-type: none">• Trimming and quality control of raw reads.• Mapping against genome reference and duplicate filter.• Variant Calling.• SNP matrix creation.• Maximum Likelihood phylogeny.• Visualization of results.
Time estimation:	1 h
Key points:	<ul style="list-style-type: none">• Importance of reference selection in SNP-based tipification.• Variant calling and SNP reconstruction is a key step in the process.• Interpretation of results is case, species and epidemiology dependant.

Introduction

Although scientific community efforts have been focused on assembly-based methods and the optimization of reconstructing complete genomes, variant calling is a essential procedure that allows per base comparison between different genomes ([Olson et al 2015](#)).

SNP-based strain typing using WGS can be performed via reference-based mapping of either reads or assembled contigs. There are many available microbial SNP pipelines, such as [Snippy](#), [NASP](#), [SNVphyl](#), [CFSAN SNP Pipeline](#), or [Lyve-SET](#).

Variant calling is a process with a bunch of potential error sources that may lead to incorrect variant calls. Identifying and resolving this incorrect calls is critical for bacterial genomics to advance. In this exercise we will use WGS-Outbreaker a SNP-based tool developed by BU-ISCI that uses bwa mapper, GATK variant caller and several SNP-filtering steps for SNP matrix contruction following maximun likelihood phylogeny using RAxML. Next image resumes the steps we are going to execute:

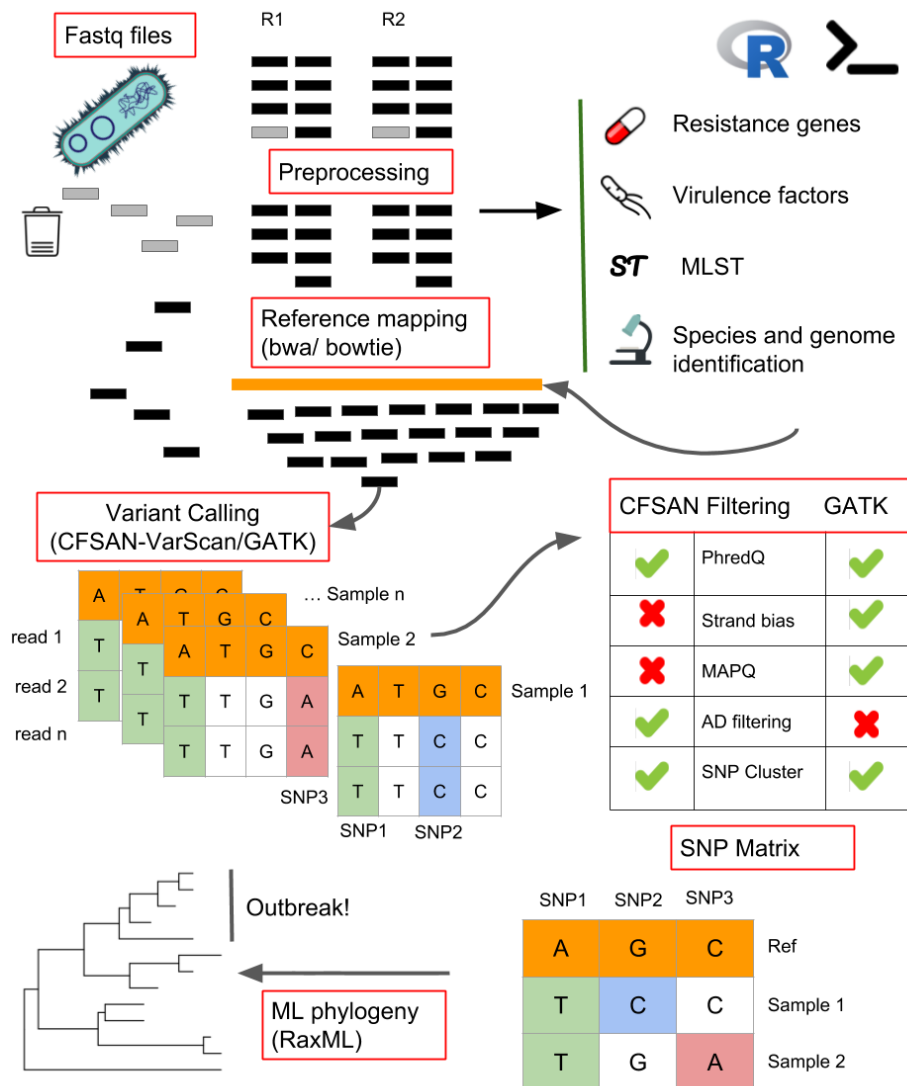

Preprocessing

We have already done our data preprocessing in the [previous exercise](#). If you remember, we executed a nextflow order which trimmed our raw reads and returned a quality report for both pre- and post-trimming files. We used FastQC for checking the data quality, Trimmomatic for the trimming and MultiQC for building the statistics report.

For mapping our reads, we will need to preprocess our data in the same way as we did for the assembly. As the results will be exactly the same we reviewed in that exercise, we will not spend more time and will move to exciting new topics.

Mapping

In the previous lecture we covered how to assemble the reads in the fastq file to recreate the original genome, or at least contigs of it. This technique requires high sequence coverage, high read length reads and good read quality, plus being highly computationally demanding. This means that it is an expensive and slow method, plus having one big disadvantage when trying to compare assembled genomes: different algorithms (and even different versions of the same software) may produce different assemblies from the same input.

For this reason, when the objective is to compare genomes of different samples, we use another method for rebuilding the genome called mapping. This technique consists in using a previously assembled genome as reference against which sequenced reads will be independently aligned against. Every read will be placed in the most likely position, ignoring any synergies between reads. This produces genomes with the same structure and coordinates that can be easily compared.

There are multiple mapping algorithms and softwares, but for this exercise we will use only [bwa](#) (H. Li and R. Durbin, 2010). bwa implements a backward search with Burrows-Wheeler Transform to efficiently align short sequencing reads against a large reference sequence such as the human genome, allowing mismatches and gaps. For longer reads, it combines its algorithm with a modified Smith-Waterman's alignment, achieving the same results as the standard algorithm but thousands of times faster. While still slower than BLAST for long query sequences, it is able to find all matches without heuristics, which makes it able to detect chimeras potentially caused by structural variations or reference misassemblies.

To map our samples with bwa, we only have to execute this command:

```
cd
cd Documents/wgs
nextflow run BU-ISCIIB/bacterial_wgs_training \
  -profile singularity \
  --reads 'training_dataset/*_R{1,2}.fastq.gz' \
  --fasta training_dataset/listeria_NC_021827.1_NoPhages.fna \
  --step mapping
```

This command will internally execute the following programs with our samples:

1. **Preprocessing** Quality control and read trimming with FastQC and Trimmomatic, as used in the previous exercise.
1. **Building bwa index** Bwa needs to build an index from the reference genome in order to now how to map the reads. This type of algorithms allows the software to do very fast searches on the genome. `bwa index fasta file`
2. **Mapping** Map each read against the reference genome using bwa mem software. `bwa mem fasta reads | samtools view -bT fasta - > bam`
3. **Post-processing and statistics** A handful of steps have to be executed before using the bam files resulting from the mapping. First, bam files have to be sorted and indexed: `samtools sort bam -o sorted.bam samtools index sorted.bam` This step allows first to reduce file size due to the compression algorithm used to generate bam files it works better if the file is sorted, and second with the file sorted and the index generated searches in the file are way faster. Imagine to look for something in a book with the pages unordered and without an index, difficult right?

Next mapping stats generated:

```
samtools stats sorted.bam > stats.txt
```

And finally we can remove some sequencing and mapping artifacts, as the duplicated reads:

```
picard MarkDuplicates \
  INPUT=sorted.bam \
  OUTPUT=dedup.bam \
  ASSUME_SORTED=true \
  REMOVE_DUPLICATES=true \
  METRICS_FILE=picardDupMetrics.txt \
  VALIDATION_STRINGENCY=LENIENT \
```

Parameters explanation:


- ASSUME_SORTED: TRUE or FALSE. Is the bam file sorted or not?
 - REMOVE_DUPLICATES: TRUE or FALSE. Do we want to remove duplicates from our bam file or do we want to just mark them?
 - METRICS_FILE: where do we want to save the metrics file?
 - VALIDATION_STRINGENCY = SILENT, LENIENT or STRICT. If an error comes out how the software is going to behave. SILENT the software just ignores the error and does not output anything, LENIENT it continues to run but outputs an error and STRICT the software stops when it encounters any error (this last option can be pretty annoying).
1. **MultiQC report** MultiQC will automatically search for the stats files and will compare them in user-friendly graphs `multiqc RESULTS_DIRECTORY` Finally, we have our mapped genomes. Now we can open them with IGV to see how they have mapped against the reference genome, and which variants are.

Mapping stats


The mapping stats are saved in `results/bwa/stats/`, where you can find a stats file per mapped sample.

We are going to visualize the summary statistics of the mapping step with MultiQC:

- The stats are represented in "Alignment metrics":

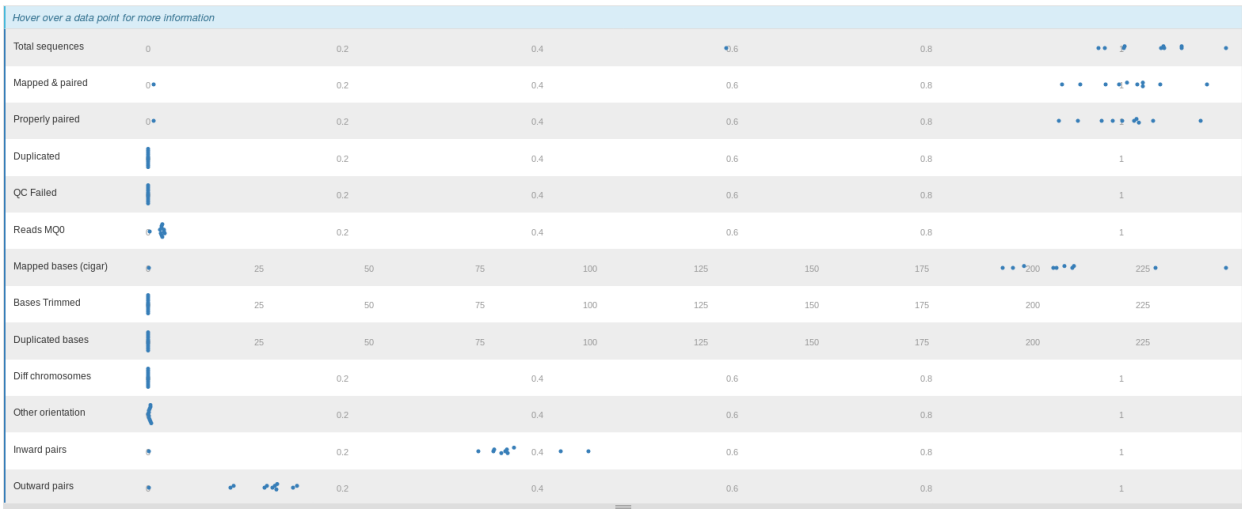


- The percentage of mapped reads per sample are plotted in "Percent Mapped":

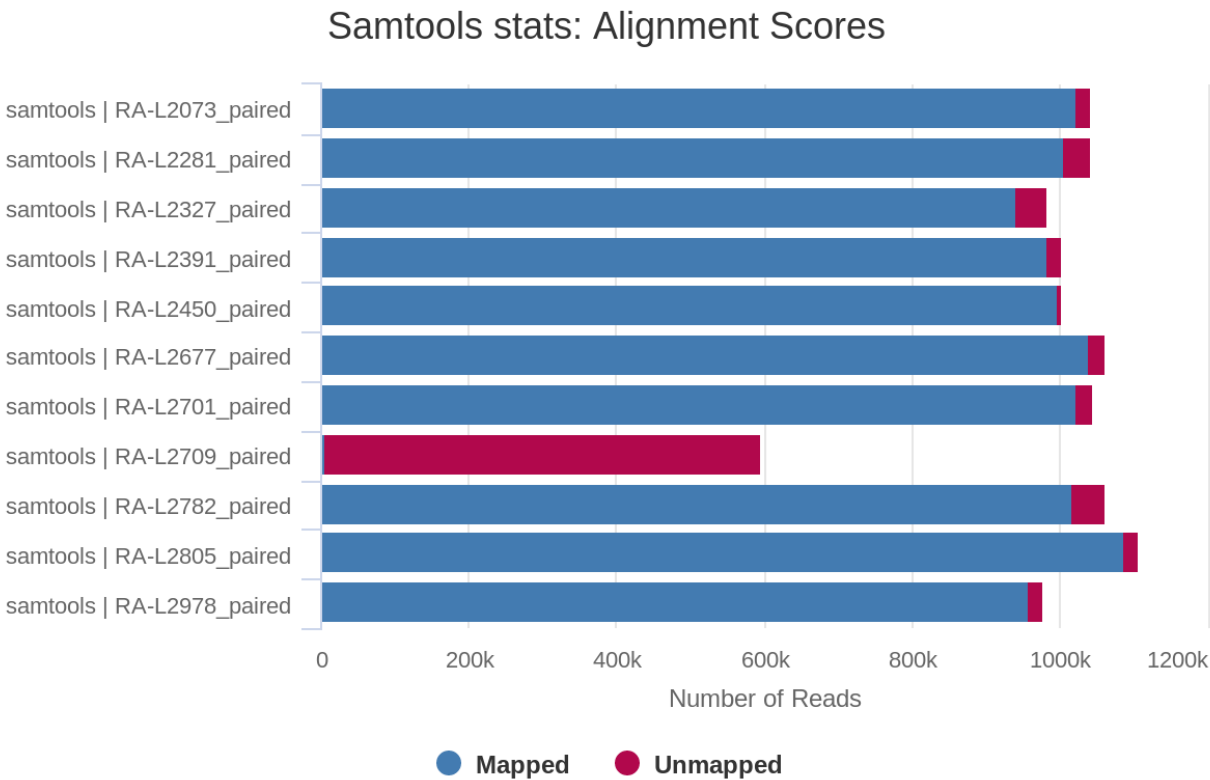

Created with MultiQC

- The number of duplicated reads per sample are plotted in "Mark Duplicates":


Picard: Deduplication Stats

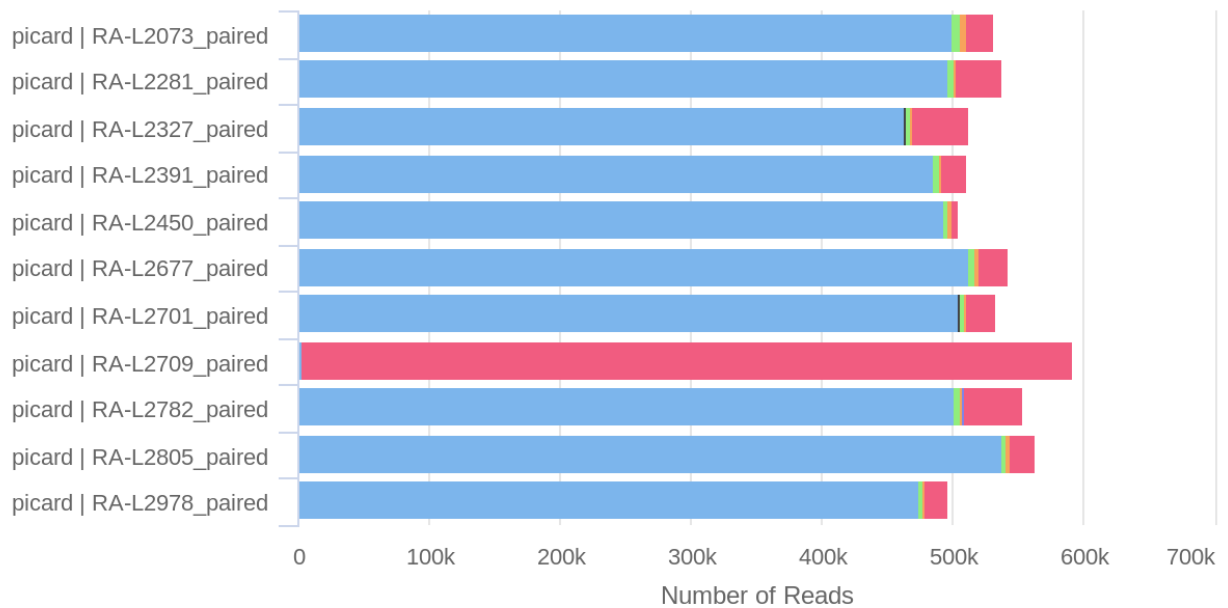

● Read Pair Unique ● Unpaired Read Unique ● Read Pair Not Optical Duplicates
● Read Pair Optical Duplicates ● Unpaired Read Duplicates ● Unmapped Reads

Created with MultiQC

Picard: Deduplication Stats

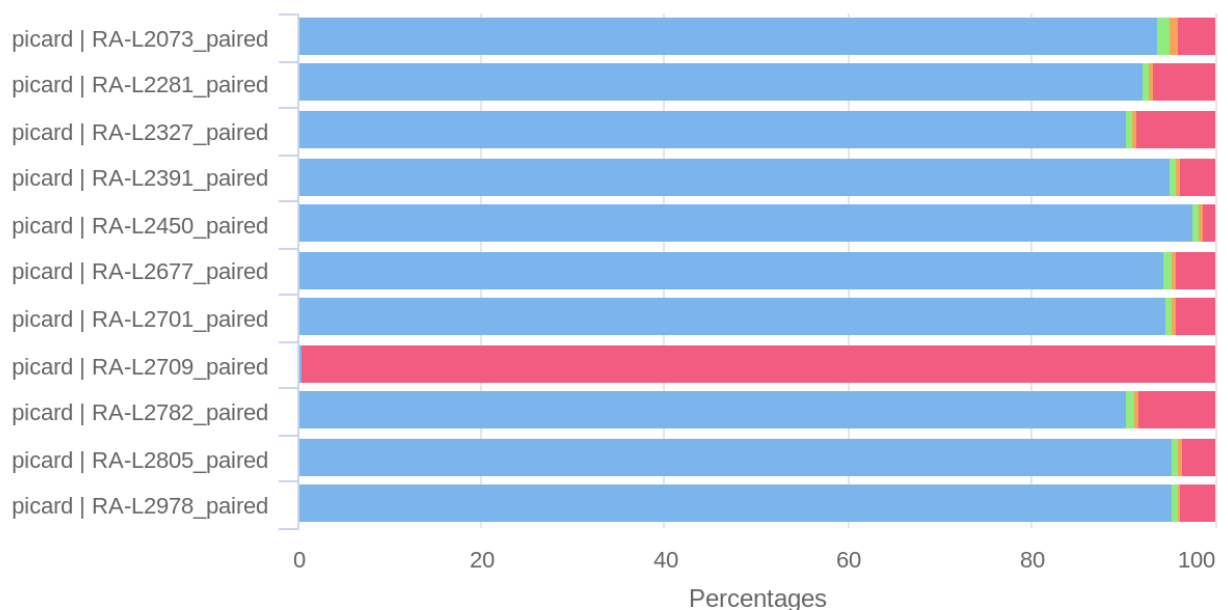

● Read Pair Unique ● Unpaired Read Unique ● Read Pair Not Optical Duplicates
● Read Pair Optical Duplicates ● Unpaired Read Duplicates ● Unmapped Reads


Created with MultiQC

Visualizing your mapping

In order to visualize our mapping we will use IGV (Integrative Genomics Viewer). This is an open source and freely available lightweight visualization tool that enables intuitive real-time exploration of diverse, large-scale genomic data sets on standard desktop computers. It supports flexible integration of a wide range of genomic data types including aligned sequence reads, mutations, copy number, RNA interference screens, gene expression, methylation and genomic annotations.


Navigation through a data set is similar to that of Google Maps, allowing the user to zoom and pan seamlessly across the genome at any level of detail from whole genome to base pair. Data sets can be loaded from local or remote sources, including cloud-based resources, enabling investigators to view their own genomic data sets alongside publicly available data.

Let's launch IGV! Navigate to your desktop and find the icon we have left for you. Double click on it and wait until it finishes loading.

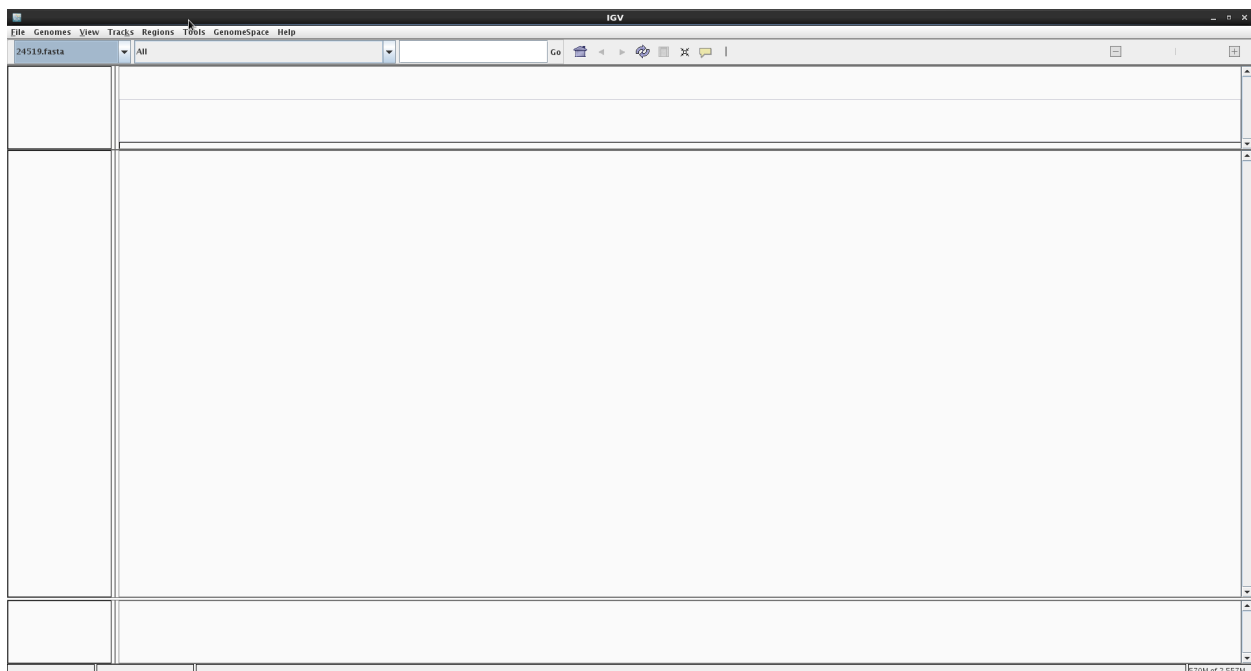

First we have to load the reference genome. Click on "Genomes" and "Load Genome from File...", navigate to our training_dataset folder and select the reference genome "listeria_NC_021827.1_NoPhages.fna".

Now, load our mapped genomes by clicking on "File" and "load from File...", navigate to our mapped genomes (HINT: they are in the shared folder inside the "results/picard" folder), and load one of them:

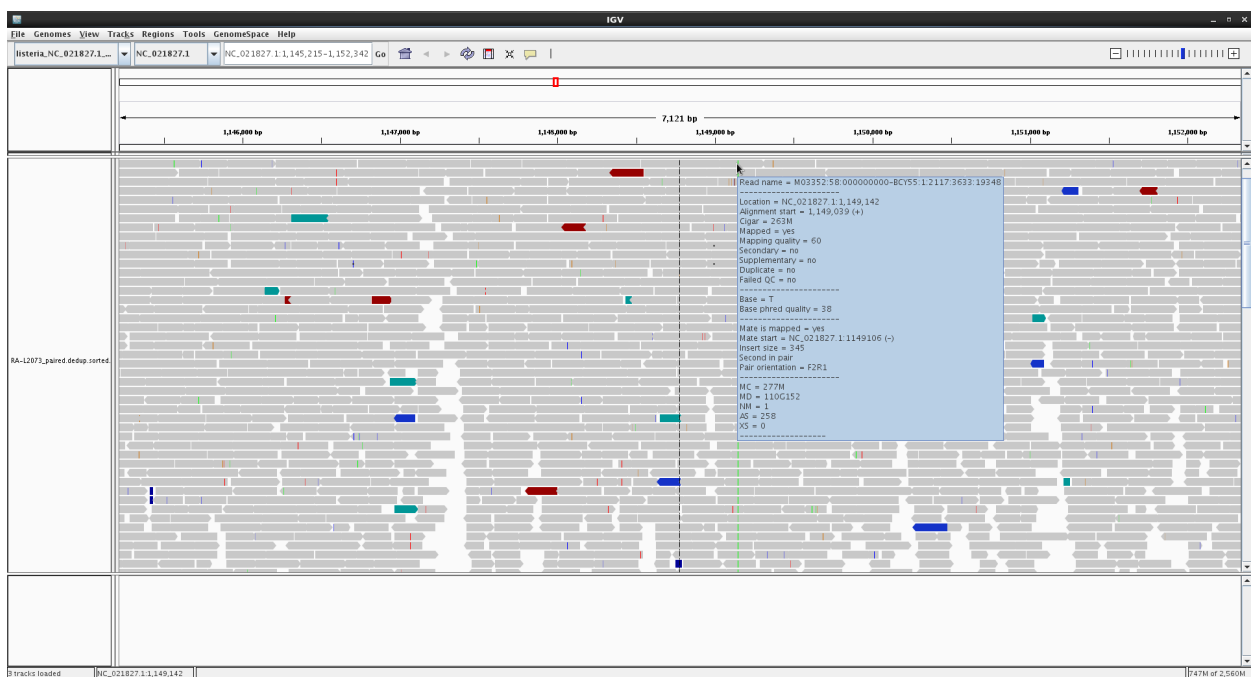

Finally, we can load as many as we want (or as many as the virtual machine survives) to compare them:



Variant Calling

We are using WGS-Outbreaker as the main software for variant calling, SNP-matrix creation and phylogeny performance. Following the development of the former exercises we are using nextflow, in this case using `outbreakSNP` step. This step includes the following processes:

- Preprocessing:
 - Trimming with trimmomatic software.
 - Quality control with fastQC.
- WGS-Outbreak software comprises the rest of steps:
 - Mapping with bwa.
 - Variant calling with GATK.
 - SNP-Matrix creation.
 - SNP-filtering:
 - PhredQ > 30
 - Strand-bias
 - MAPQ
 - SNP cluster, < 3 SNPs / 1000 pb

Everything clear..? So let's run it.

Run the exercise

First of all we need to be clear in which folder we are. We need to be in our working directory `/home/alumno/Documents/wgs` and our training dataset downloaded the first day must be there (If you had any problem the previous sessions please refer to the [setup tutorial](#)).

You can run this command to check where you are:

```
pwd
```

Output:

```
/home/alumno/Documents/wgs
```

And this one to list all the files in your working directory. Check there is the training_dataset folder and the results folder from previous sessions.

```
ls
```

Output:

```
training_dataset results work
```

Once our localization is correct we will launch nextflow with the next parameters:

- Raw reads
- step outbreakSNP
- saveTrimmed -> this parameters saves the fastq files trimmed in our results dir.

- `outbreaker_config` <- config file with all the parameters required by WGS-Outbreaker

```
nextflow run BU-ISCIIII/bacterial_wgs_training \
  -profile singularity \
  --reads 'training_dataset/*_R{1,2}.fastq.gz' \
  --fasta training_dataset/listeria_NC_021827.1_NoPhagues.fna \
  --step outbreakSNP \
  --saveTrimmed \
  --outbreaker_config training_dataset/wgsoutbreaker.config
```

Output:

```
N E X T F L O W ~ version 0.29.0
Launching `main.nf` [distracted_magritte] - revision: 3508cbd2da
WARN: Process `multiqc` is defined two or more times
WARN: Process `multiqc` is defined two or more times
WARN: Process `multiqc` is defined two or more times
=====
BU-ISCIIII/bacterial_wgs_training : WGS analysis practice v1.0
=====
Reads                : test/full_dataset/*_R{1,2}*.fastq.gz
Data Type            : Paired-End
Fasta Ref            : test/listeria_NC_021827.1_NoPhagues.fna
Keep Duplicates      : false
Step                 : outbreakSNP
Container            : ./wgs_bacterial.simg
Current home         : /home/smonzon
Current user         : smonzon
Current path         : /home/smonzon/Documents/desarrollo/bacterial_wgs_training
Working dir          : /home/smonzon/Documents/desarrollo/bacterial_wgs_training/work
Output dir           : results
Script dir           : /home/smonzon/Documents/desarrollo/bacterial_wgs_training
Save Reference       : false
Save Trimmed         : true
Save Intermeds       : false
Trimmomatic adapters file: $TRIMMOMATIC_PATH/adapters/NexteraPE-PE.fa
Trimmomatic adapters parameters: 2:30:10
Trimmomatic window length: 4
Trimmomatic window value: 20
Trimmomatic minimum length: 50
Config Profile       : singularity
=====
[warm up] executor > local
[ca/cbb117] Submitted process > fastqc (RA-L2805)
[4b/65f7a1] Submitted process > fastqc (RA-L2450)
[32/8ebe88] Submitted process > fastqc (RA-L2281)
[23/04ab41] Submitted process > fastqc (RA-L2073)
[a7/f9e938] Submitted process > fastqc (RA-L2391)
[75/709471] Submitted process > trimming (RA-L2073)
[94/87b2b5] Submitted process > makeBWAindex (listeria_NC_021827.1_NoPhagues)
[94/b39b86] Submitted process > trimming (RA-L2805)
[df/e01505] Submitted process > fastqc (RA-L2709)
.....
BU-ISCIIII - Pipeline complete
```

This will take a while so we need to move forward and understand what we are doing and learn how to see and interpret our results.

Understanding WGS-Outbreaker config file

First of all, let's take a look to the config file for a moment: [WGS-Outbreaker config file](#). This file will allow us to configure all necessary parameters for running WGS-Outbreaker. The file is organized in several sections. 1. **Steps configuration**: in this section we can select with YES/NO which pipeline steps we would want to run, in this case we have prefilled the steps that we can run in this training.


```
##### Pipeline steps: Fill in with YES or NO (capital letter) #####
TRIMMING=NO
CHECK_REFERENCES=YES
MAPPING=YES
DUPLICATE_FILTER=YES
VARIANT_CALLING=YES
KMERFINDER=NO
SRST2=NO
CFSAN=NO
VCF_TO_MSA=YES
RAXML=YES
STATS=YES
```

1. **Input data:** we can provide the path where the input files are, and the path where we want our results. Also we can include our sample names and the raw reads filenames. These reads must be in the input directory provided. `` # Directory with input files
INPUT_DIR=/home/smonzon/Documents/desarrollo/bacterial_wgs_training/results/trimming

Directory for output files

OUTPUT_DIR=/home/smonzon/Documents/desarrollo/bacterial_wgs_training/results/wgs_outbreaker

```
##### INPUT VARIABLES#####
```

Samples info:

All samples ID must be separated by ":", then for each sample there must be a line with the names for

R1 and R2 separated by tabulator

Example:

```
#=AAAA_01:BBBB_02
# AAAA_01=AAAA_01_R1.fastq.gz      AAAA_01_R2.fastq.gz
# BBBB_02=BBBB_02_R1.fastq.gz      BBBB_02_R2.fastq.gz
```

Moreover we have to include the path for our reference files:

```
##### Reference Variables #####
```

Path to reference genome

GENOME_REF=listeria_NC_021827.1_NoPhages.fna

Path to reference genome without ".fasta"

GENOME_NAME=listeria_NC_021827.1_NoPhages

```
3. **Trimming, mapping, variant calling and phylogeny parameters:** The end of the config file includes a series of default values. For example the most variable parameter we can probably find is the maximum number of SNPs we are going to allow in a window.
```

```
##### SNP FILTERS #####
```

The maximum number of SNPs allowed in a window.

MAX_SNP=3

The length of the window in which the number of SNPs should be no more than max_num_snp

WINDOW_SIZE=1000

Results analysis

Let's proceed to analyze the results. We can find them in:

/home/alumno/course_shared_folder/wgs_outbreaker

This directory contains several folders including:

└─ Alignment -> already analyzed └─ QC -> already analyzed └─ raw -> symbolic links to raw reads └─ RAXML -> phylogenetic results
└─ stats -> alignment and variant calling stats. └─ variant_calling -> variant calling files.

Since alignment and quality control results has been previously addressed in this course (see [02_QualityAndAssembly.md]) (

Variant calling results

Variants are stored in plain text files in vcf format (variant calling format). Vcf files can be found in:

wgs_outbreaker/variant_calling/variants_gatk/variants

Here we can find a bunch of vcf files for each filtering steps we made:

- *.g.vcf <- this file contains a special vcf format that includes both variant and invariants sites information.
- snps_indels.vcf <- contains raw variants, both indels and snps found by GATK in the samples. This is a multisample vcf
- In order to follow GATK's best practice protocol for high quality variant filtering, snps and indels must be treated s

```
##fileformat=VCFv4.2 ##FILTER= ##FILTER= ##FILTER= ##FILTER= ##FILTER= ##FILTER= 3.0"> ##FORMAT= ##FORMAT= ##FORMAT=
##FORMAT= ##INFO= ##INFO= ##INFO= ##INFO= ##INFO= ##INFO= ##INFO= ##INFO= ##contig= ##contig=
##reference=file:///home/smonzon/Documents/bacterial_wgs_training/work/ba/20837f0b9838403205e62589d7ac8b/listeria_NC_021827.1_NoPhagues.fna
##source=SelectVariants #CHROM POS ID REF ALT QUAL FILTER INFO FORMAT RA-L2073 NC_021827.1 276 . C A 291.68
RMSMappingQuality;SnpCluster;StandOddRatio
AC=1;AF=0.100;AN=10;DP=35;FS=0.000;MLEAC=1;MLEAF=0.100;MQ=38.49;QD=34.24;SOR=4.407 GT:AD:DP:GQ:PL 0:3,0:3:99:0,117
NC_021827.1 731 . A G 2313.68 SnpCluster
AC=1;AF=0.100;AN=10;DP=101;FS=0.000;MLEAC=1;MLEAF=0.100;MQ=60.00;QD=33.05;SOR=0.811 GT:AD:DP:GQ:PL 0:3,0:3:99:0,101
NC_021827.1 921 . C T 1841.68 SnpCluster
AC=1;AF=0.100;AN=10;DP=110;FS=0.000;MLEAC=1;MLEAF=0.100;MQ=60.00;QD=29.70;SOR=0.826 GT:AD:DP:GQ:PL 0:3,0:3:99:0,101
NC_021827.1 1067 . C T 2250.68 SnpCluster
AC=1;AF=0.100;AN=10;DP=168;FS=0.000;MLEAC=1;MLEAF=0.100;MQ=60.00;QD=31.70;SOR=0.779 GT:AD:DP:GQ:PL 0:3,0:3:99:0,101
NC_021827.1 2114 . C T 8324.89 SnpCluster
AC=3;AF=0.300;AN=10;DP=341;FS=0.000;MLEAC=3;MLEAF=0.300;MQ=60.00;QD=30.49;SOR=0.841 GT:AD:DP:GQ:PL 0:3,0:3:99:0,101
NC_021827.1 2180 . G A 7855.89 SnpCluster
AC=3;AF=0.300;AN=10;DP=342;FS=0.000;MLEAC=3;MLEAF=0.300;MQ=60.00;QD=28.67;SOR=0.832 GT:AD:DP:GQ:PL 0:3,0:3:99:0,101
```

- Finally we continue to filter snps calls for our SNP matrix, and we filter SNPs which are included in a window of 1000

Phylogeny results

Phylogenetic tree reconstruction is performed using RAXML with 100 inferences and 100 bootstrap repetitions. RAXML resul

/home/Alumno/Documents/wgs/results/RAXML/{variant_caller}

Two different trees are generated one with only SNPs passing all filters (preser) and one with all snps (all_snp). Both

RAXML outputs one file per inference and per bootstrap so the folder is full of files. Don't worry we only need the fina

RAXML_bipartitions.RAXML_TREE_ANNOT

Now we are going to open firefox browser and go to [iTOL website](https://itol.embl.de/). This web allows us to visualiz

Once in iTOL website we click in Upload in the top menu. Next we upload our tree as shown in the image:

<p align="center"></p>

Now we are visualizing our tree and we can manipulate it. First of all, as we are facing an unrooted tree with long and

<p align="center"></p>

iTOL offers multiple annotation and manipulation options. We can selecto for example the display of bootstrap in the adv

<p align="center"></p>

You can play with all the options iTOL offers reading its [documentation](https://itol.embl.de/help.cgi), and export the

Now we are going to focus on our results. Our final tree should look something like this:

<p align="center"></p>

Which strains do you think belong to the outbreak?

Following the instructions seen in the theory class, here we can focus on topology and bootstrap.

Do we find any monophyletic group?

Does it have more than 80 bootstrap value?

SNPs distance

As we have studied in the theory class, maximum likelihood methods for phylogeny only offers as branch lenght the averag

In order to know exactly the SNPs differences between the strains WGS-Outbreaker outputs a distance matrix showing the p

We can check this file in:



/home/alumno/Documents/wgs/results/wgs_outbreaker/stats

```
|dis_matrix.names|RA.L2073|RA.L2281|RA.L2327|RA.L2391|RA.L2450|RA.L2677|RA.L2701|RA.L2782|RA.L2805|RA.L2978|
|-----|-----|-----|-----|-----|-----|-----|-----|-----|-----|
|RA-L2073|0|9403|9028|80|46|46|49|9120|4|2|
|RA-L2281|9403|0|8777|9415|9397|9397|9402|9183|9403|9401|
|RA-L2327|9028|8777|0|9040|9022|9022|9027|4277|9028|9026|
|RA-L2391|80|9415|9040|0|74|74|79|9132|80|78|
|RA-L2450|46|9397|9022|74|0|38|45|9114|46|44|
|RA-L2677|46|9397|9022|74|38|0|45|9114|46|44|
|RA-L2701|49|9402|9027|79|45|45|0|9119|49|47|
|RA-L2782|9120|9183|4277|9132|9114|9114|9119|0|9120|9118|
|RA-L2805|4|9403|9028|80|46|46|49|9120|0|2|
|RA-L2978|2|9401|9026|78|44|44|47|9118|2|0|
```

As we see the SNP difference cutoff is important here, and it will depend on the strain and the case. If we stablish 3-5

